


Fermentation Monitoring

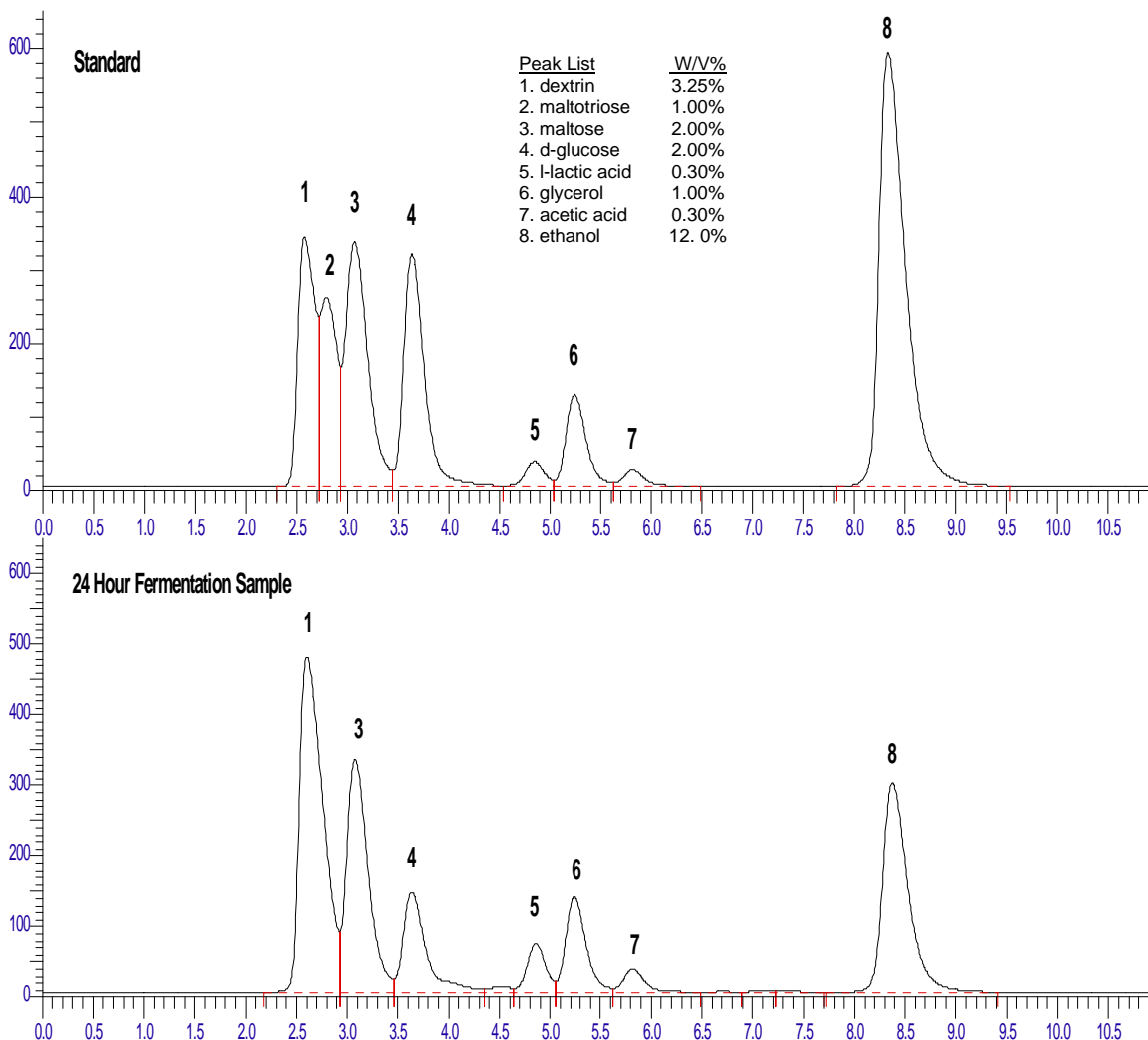

Conditions

Mobile phase: 0.001M H₂SO₄
 Flow: 0.8 mL/min
 Temperature: 60°C
 Detector: Refractive index @ 40°C
 Injection volume: 10 µL

Column

BIO-RAD column: Aminex® Fermentation Monitor
 Dimensions: 150 x 7.8 mm
 Part number: 1250115

PerkinElmer Life and Analytical Sciences

710 Bridgeport Avenue
 Shelton, CT 06484-4794 USA
 Phone: (800) 762-4000 or
 (+1) 203-925-4602
www.perkinelmer.com

