


Research Use Only. Not for use in diagnostic procedures.

## Cysteine, L-[<sup>35</sup>S]-

Product Number: NEG022T

### LOT SPECIFIC INFORMATION

Lot Number:	_____01130_____
Specific Activity:	_____1075.0_____ Ci/mmol
	_____39.8_____ TBq/mmol
Concentration:	_____11.0_____ mCi/ml
	_____407.0_____ MBq/ml
	_____0.0102_____ μmol/ml
	_____0.0012_____ mg/ml
Calibration Date:	_____7-Feb-2020_____

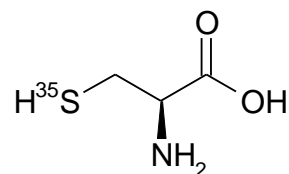

M.W. 121.15

C<sub>3</sub>H<sub>7</sub>NO<sub>2</sub>S

**PACKAGING:** Aqueous solution containing 10mM dithiothreitol in a NENSURE vial with polypropylene-v-insert, shipped in dry ice.

### STABILITY AND STORAGE:

- The rate of decomposition of Cysteine, L-[<sup>35</sup>S]- is less than 3% per month when stored as supplied at -80°C. Storage at higher temperatures will result accelerated breakdown.
- Regardless of storage temperature, repeated freezing and thawing can result in product degradation. If the entire quantity will not be used at once, it is recommended that the product be thawed quickly, aliquoted under N<sub>2</sub> into samples of an appropriate size, refrozen immediately and stored at -20°C or below.
- Removal of DTT is not advised. This results in rapid oxidation of the material.

**HAZARD INFORMATION:** WARNING: This product contains a chemical known to the state of California to cause cancer.

### QUALITY CONTROL:

**Radiochemical Purity:** This lot was initially found to be >95% when determined by high pressure liquid chromatography on a C18 column.

**Chemical Purity:** As determined by amino acid analysis, NEG-022T contains the unlabeled amino acid serine. This may be present in up to a 3:1 molar ratio relative to Cysteine, L-[<sup>35</sup>S].

**PREPARATIVE PROCEDURE:** Cysteine, L-[<sup>35</sup>S] is obtained from the protein hydrolysate of bacteria grown in the presence of carrier-free [<sup>35</sup>S] sulfate. The product is purified by high pressure liquid chromatography.

**SAFE HANDLING:** Because this product has been stored at -80°C, it is possible that pressure may develop in the vial during the thawing process. In addition, volatile <sup>35</sup>S-labeled decomposition products are generated at a rate of ~0.01% per week. We recommend that prior to opening, vials are first vented in a fume hood using the following procedure:

- Slide aside the dust cover on the cap to expose the septum.
- Pierce the septum with a cotton-plugged syringe needle or charcoal trap (NENTM NEX-033T), taking care that the tip does not come in contact with the product.
- If the product is frozen, quickly thaw at room temperature or in a 37°C water bath. Any pressure developed will vent through the syringe needle.
- Remove the needle and dispose of as contaminated equipment.

When used for *in vivo* labeling experiments, we recommend that specific steps be taken to minimize incubator and water bath contamination. We suggest using a shallow tray of activated charcoal, charcoal sticks or charcoal filter units to trap <sup>35</sup>S volatiles and reduce contamination.

### SPECIAL INFORMATION:

Decay of <sup>35</sup>S (physical half-life, 87.4 days):

		DECAY FACTORS									
Days BEFORE	Assay Date	0	1	2	3	4	5	6	7	8	9
	30	1.269	1.279	1.289	1.299	1.309	1.320	1.330	1.341	1.352	1.362
	20	1.172	1.181	1.191	1.200	1.210	1.219	1.229	1.239	1.249	1.259
	10	1.083	1.091	1.100	1.109	1.117	1.135	1.135	1.144	1.153	1.163
	0	1.000	1.008	1.016	1.024	1.032	1.049	1.049	1.057	1.066	1.074
Days AFTER	Assay Date	0	1	2	3	4	5	6	7	8	9
	0	1.000	0.992	0.984	0.976	0.969	0.961	0.954	0.946	0.939	0.931
	10	0.924	0.916	0.909	0.902	0.895	0.888	0.881	0.874	0.867	0.860
	20	0.853	0.847	0.840	0.833	0.827	0.820	0.814	0.807	0.801	0.795
	30	0.788	0.782	0.776	0.770	0.764	0.758	0.752	0.746	0.740	0.734

The specific activity at time t (SA<sub>t</sub>) may be calculated, using the following equation, from the specific activity at the calibration date (SA<sub>0</sub>) and the decay factor (f) given above.

$$SA_t = \frac{f}{1/SA_0 - (1-f)/1494}$$

**REFERENCE:** Rubin, I.B., and Goldstein, G., (1970) *Anal. Biochem.* 33, 244-254.

### RELATED PRODUCTS:

NEG-009A      Methionine, L-[<sup>35</sup>S]-  
NEG-009C      Methionine, L-[<sup>35</sup>S]-  
NEG-009H      Methionine, L-[<sup>35</sup>S]-  
NEG-009T      Methionine, L-[<sup>35</sup>S]-  
1084:121007

PerkinElmer, Inc.  
549 Albany Street  
Boston, MA 02118 USA  
P: (800) 762-4000 or (+1) 203-925-4602  
[www.perkinelmer.com/enradiochemicals](http://www.perkinelmer.com/enradiochemicals)

For a complete listing of our global offices, visit [www.perkinelmer.com/ContactUs](http://www.perkinelmer.com/ContactUs)  
Copyright ©2010, PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc.  
All other trademarks are the property of their respective owners.

