

Caution: For Laboratory Use. A research chemical for research purposes only

Eu-W1024-labeled Anti-Phospho-eIF4E-binding protein 1 (Thr37/46) Antibody

Product No.: TRF0216-D / TRF0216-M

Lot No.: 2854433

Material Provided

Format: TRF0216-D 10 µg (1 562 assay points*)
TRF0216-M 100 µg (15 625 assay points*)
*Assuming 40 fmol/assay point

Manufacturing Date: April 8, 2021

Product Information

Antibody: Europium-labeled rabbit monoclonal antibody recognizing human eukaryotic translation initiation factor 4E-binding protein 1 phosphorylated at Thr37 or Thr46.

Storage Buffer: 50 mM Tris-HCl (pH 7.4), 0.9% NaCl with 0.05% sodium azide as preservative and 0.1 % BSA

Molecular Weight: 160 000

Stability: This product is stable for at least **18 months** from the manufacturing date when stored in its original packaging and at the recommended storage conditions.

Storage Conditions: Store at 4°C.

Safety Note: The storage buffer contains sodium azide (NaN₃) as a preservative. Disposal of all waste should be in accordance with local regulations.

Quality Control

This product passed QC release specifications based on spectrophotometric analysis of the labeled antibody.

Labeling Ratio: 6.79 (Eu/Ab)

Concentration: 100 µg/mL (0.625 µM)

Recommended Assay Conditions

mTOR KINASE: ATP TITRATION

SUGGESTED METHOD:

(Specific applications might require optimization)

Reagent Preparation

- Prepare 1X Kinase Assay Buffer: 50 mM HEPES pH 7.5, 1 mM EGTA, 10 mM MgCl₂, 3 mM MnCl₂, 2 mM DTT and 0.01% Tween-20.
- Prepare a 4X *ULight*-eIF4E-binding protein 1 (Thr37/46) Peptide solution: dilute *ULight*-eIF4E-binding protein 1 to a concentration of 200 nM in Kinase Assay Buffer.
- Prepare a 2X mTOR solution: dilute enzyme to a concentration of 1 nM in Kinase Assay Buffer. Keep on ice.
- Prepare a 4X ATP solution: dilute ATP to concentrations ranging from 40 nM to 4 mM (serial half-log dilutions) in Kinase Assay Buffer. Keep on ice.
- Prepare 1X Detection Buffer: dilute 0.5 mL of 10X Detection Buffer with 4.5 mL of H₂O.
- Prepare a 4X Stop solution*: dilute EDTA to a concentration of 32 mM in 1X Detection Buffer.
- Prepare a 4X Detection Mix: dilute Europium-anti-phospho-eIF4E-binding protein 1 (Thr37/46) Antibody to a concentration of 8 nM in 1X Detection Buffer.


Protocol

- Pipet 5 µL of 2X mTOR solution (0.5 nM final concentration).
- Add 2.5 µL of 4X *ULight*-eIF4E-binding protein 1 (Thr37/46) peptide solution into a 384-well white OptiPlate-384 (50 nM final concentration).
- Add 2.5 µL of 4X ATP solution (10 nM to 1 mM final concentrations).
- Cover plate with TopSeal-A and incubate 60 min at 23°C.
- Add 5 µL of 4X Stop solution* and incubate 5 min at 23°C.
- Add 5 µL of Detection Mix (2 nM Eu-anti-phospho-eIF4E-binding protein 1 (Thr37/46) Antibody final concentration).
- Cover plate with TopSeal-A and incubate 60 min at 23°C.
- Remove TopSeal-A and read in TR-FRET mode at 665 nm (excitation at 320 or 340 nm).

**Alternatively, the Stop solution and Detection Mix can be premixed and added together to the kinase reaction.*

TYPICAL ATP TITRATION DATA obtained using the EnVision® Multilabel Reader:

mTOR kinase assay using *ULight*[™]-4E-BP1 (Thr37/46) Peptide and Eu-anti-phospho-4E-BP1 (Thr37/46) Antibody

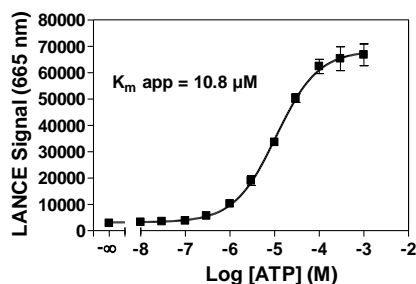

Suggested Materials

- Substrate: *ULight*[™]-4E-BP1 (Thr37/46) Peptide
- Antibody: Europium-anti-P-4E-BP1 (Thr37/46)
- Kinase: mTOR, active
- Detection Buffer: LANCE[®] Detection Buffer, 10X
- Plate: OptiPlate[™]-384, white
- TopSeal[™]-A

Supplier

Cat. No.

PerkinElmer	TRF0128
PerkinElmer	TRF0216
Millipore/Upstate	14-770
PerkinElmer	CR97-100
PerkinElmer	6007299
PerkinElmer	6050195

Please visit our website for additional resource:

www.perkinelmer.com/LANCE

This product is not for resale or distribution except by authorized distributors.

LANCE[®] products may be covered by one or more the following patents: US 5,202,423 and US 5,324,825

LIMITED WARRANTY: PerkinElmer BioSignal Inc. warrants that, at the time of shipment, the products sold by it are free from defects in material and workmanship and conform to specifications which accompany the product. PerkinElmer BioSignal Inc. makes no other warranty, express or implied with respect to the products, including any warranty of merchantability or fitness for any particular purpose. Notification of any breach of warranty must be made within 60 days of receipt unless otherwise provided in writing by PerkinElmer BioSignal Inc. No claim shall be honored if the customer fails to notify PerkinElmer BioSignal Inc. within the period specified. The sole and exclusive remedy of the customer for any liability of PerkinElmer BioSignal Inc. of any kind including liability based upon warranty (express or implied whether contained herein or elsewhere), strict liability contract or otherwise is limited to the replacement of the goods or the refunds of the invoice price of goods. PerkinElmer BioSignal Inc. shall not in any case be liable for special, incidental or consequential damages of any kind.

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com



For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2009, PerkinElmer, Inc. All rights reserved. PerkinElmer[®] is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.